

Instructions

The following instructions will allow you to receive SMS Text Messages on your phone from the Comark Cloud.

Setup a Text Magic Account

If you already have a Text Magic Account you can skip this section.

Visit www.textmagic.com and register for an Account

Log into TextMagic and Click on Services – Distribution Lists to bring up this page.....

Distribution lists

[Watch tutorial](#)

All emails sent to your SMS distribution list addresses are immediately forwarded to the mobile phones of distribution list members as SMS.

+ New distribution list	Delete	<input type="text" value="Search"/>	<input type="button" value="Q"/>	
<input type="checkbox"/>	Name ▾	Email for sharing ▾	Recipients ▾	Actions

Click on New distribution list.

Complete the information on the page....

New distribution list

1. Add details **2. Finish up**

Distribution list name

Add recipients x
Start typing a number or contact name

[Contacts](#)
[Lists](#)
[Frequently sent](#)

Recipients: 1

Include in SMS content

Email subject and email body
The email subject and body will be included in SMS content.

Only email body
Only the email body will be included in SMS content.

Only email subject
Only the email subject will be included in SMS content.

Maximum length

Send from

[Cancel](#) **Continue** ▶

Click on Continue to finish.


You will get an e-mail address which you need to copy.

Add a new user to the Cloud


You now need to add a new user to the Comark Cloud and associate that user with the Locations/Loggers that you wish to send SMS alerts from.

Log on to the Cloud (You must be an Account Administrator in order to Add more Users to the Cloud)

Click on Account (top right of the screen) and then Click on the USERS Tab....

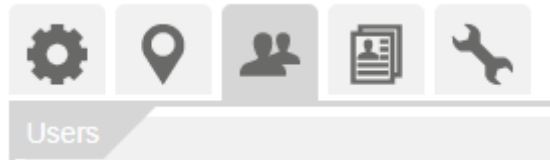

At the bottom of the page, select ADD and you will be presented with this page....

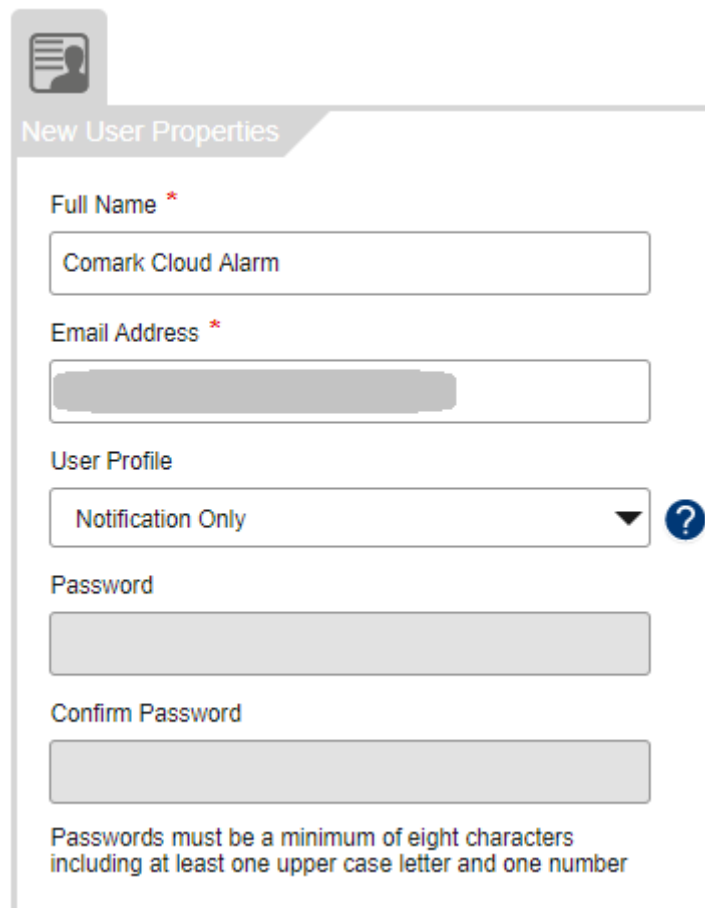A screenshot of a web form titled "New User Properties". It contains several input fields: "Full Name *" with the value "Comark Cloud Alarm"; "Email Address *" which is blurred; "User Profile" with a dropdown menu set to "Notification Only" and a help icon; "Password" and "Confirm Password" fields, both blurred. At the bottom, there is a note: "Passwords must be a minimum of eight characters including at least one upper case letter and one number".

Select Notification Only for the User.

Now Select the Locations for the Alarms SMSs.

Click on Apply to Save settings.

Allocate Transmitters to Receive SMS Alerts


Each Transmitter on the Cloud that you want to send SMS Alerts has to be setup.

From the Home Page, select Devices. Navigate to the device you want to setup.




Tick the box and then select the

Then select the Notifications Tab, see image below.

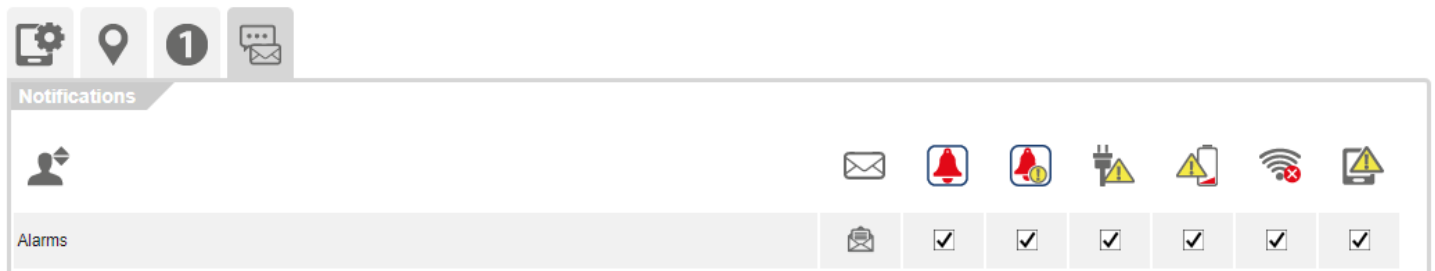

From here simply select the User as previously setup and then tick the appropriate boxes depending on the type of Alerts you wish to receive as SMS.

Note: The type of E-mail for SMS must be in Text Format, see picture below.

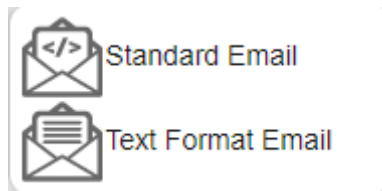

The setup process is finished.

End.